

# **MEDIT** *i* 700 wireless



# Magic Made Easy with a Simple Touch

Get the power to take your practice to the next level!



# Beyond i700's Proven Tech

Still the same **super fast**, **light**, and **accurate** *i***700**. But now **wireless**.



### **No More Wires**

Free from wires, scanning is comfortable at any angle. Scan proximal areas effortlessly without constraints in movement.



#### Fast, Accurate. Connected

Scanning speed of up to 70 FPS\*. You can experience the amazingly smooth and fast scanning performance of Medit's best-selling model. Enjoy the fast speed of the *i*700 wirelessly, flawlessly.



#### **Supercharged Yet Light**

Outstanding weight balance even with a battery attached for added comfort.



## Powerful Hardware Meets Versatile Software

Optimized for convenience. Not only from hardware aspects, but also software. *i*700 wireless, packed with the latest of MEDIT technology.

# Full arch accuracy: 10.9 µm ± 0.98 Up to 70 FPS

## **Reversible Tip**

Small Tip (Sold separately)

For easy application to the posterior region for patients with small oral cavity No need to twist your wrist for maxillary scans

### Wider Scan Area

Enlarged scanning window captures large areas at once while maintaining the small size of the scanner head



the handpiece

## **Auto Calibration**

Automatic tool makes calibration a breeze



Auto Calibration Tool (Sold separately) When the scanner's cable disappears, so do restrictions to scan movements. Strengthened impact resistance minimizes damage from accidental drops. Weight balanced design of the *i*700 wireless provides a stable, comfortable scanning experience.




# The Best Battery Life Ever

## 1 hour of continuous scanning, up to 8 hours in standby mode



All-Day-Scan-with-Single-Battery. Scan all day without worrying about charging.

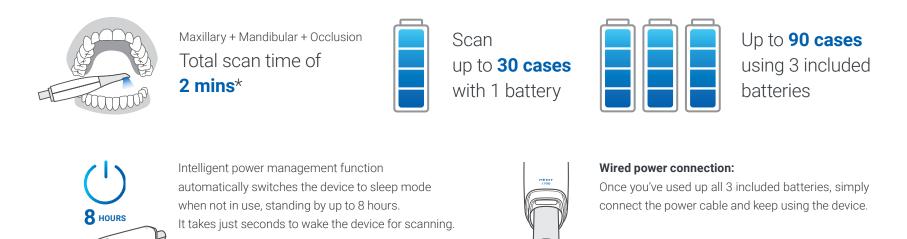

# Magical Scanning, Easier Than Ever



## Experience amazing digital dentistry with **Medit Scan for Clinics**. Our cutting-edge software technology makes scanning more convenient for you.

## For scanning assistance



help for abnormalities detected during eases the scanning process by filtering scanning.

out soft tissue such as lips, cheeks, and tongue.

Automatically filters registered colors, such as the color of medical gloves, during scanning.

Scan with your own scanning strategy. Fragmented scan data are aligned with automatically detected common area.

## For more precise and accurate scan data

### Smart Shade Guide



Automatically distinguishes and suggests the best matching shade from scan data.

A.I. Scan Body Matching 16 💐 🐒 就 to 🎘 🔹 🕱 🔍 8

Secures accurate scan body data in Secures accurate margin data in difficult-to-scan areas, such as the difficult-to-scan situations. posterior region.

A.I. Abutment Matching N 10 📰 

HD Mode




**HD** HD Off On O

Get a sharper scan image in the margin area with the upgraded HD Mode.

# Medit Link

Begin your clinic's data management with Medit Link.

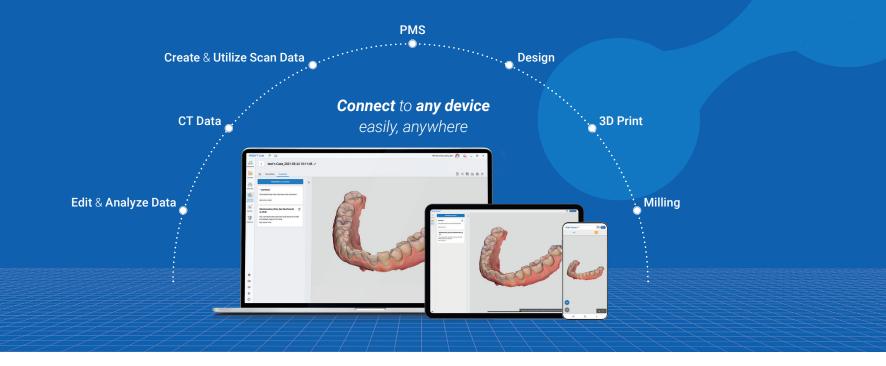

You can create and manage scan data matched with patient information in the electronic chart program and link them with treatment solutions of various workflows. The various apps of the Medit Link App Box help provide higher standard consultation that is easier to understand for patients.

Medit Link's cloud storage encrypts and stores data for reliable protection against virus and ransomware attacks. It also enables data sharing and real-time communication between dentists, patients, and lab technicians anytime, anywhere on a PC, a tablet, or any mobile device with an internet connection.

## Digital Dentistry with Medit Link & Medit Apps

Digitize Your Treatment with Medit Link



## **Use Your Scan Data for Much More**

Scan data created with the Medit *i*700 wireless can be linked to the workflow of various solution providers via Medit Link. The data can also be used with any software that can open standard 3D formats such as STL, PLY, and OBJ. And the wide utilization of our open system is expanding as we speak.


#### What you can do in Medit Link :

Clear aligners | Crown and bridges | Implant-supported restorations | Full and partial dentures | Orthodontics | Implant surgical guides | Sleep apnea devices | Dental models | 3D printing...



## Medit Apps

Medit Apps quickly introduce digital workflows using scan data and help you to utilize them in a diverse range of cases. You can analyze data, create restorations and models, and communicate with your patients in the simplest way. With Medit Apps, a lot more can be done on the chairside like designing splints, temporaries, bridges, and so on! By the way, you don't need to have a Medit scanner or pay a huge annual fee to join the Medit Ecosystem. Feel literally free to enjoy digital dentistry with Medit Apps.

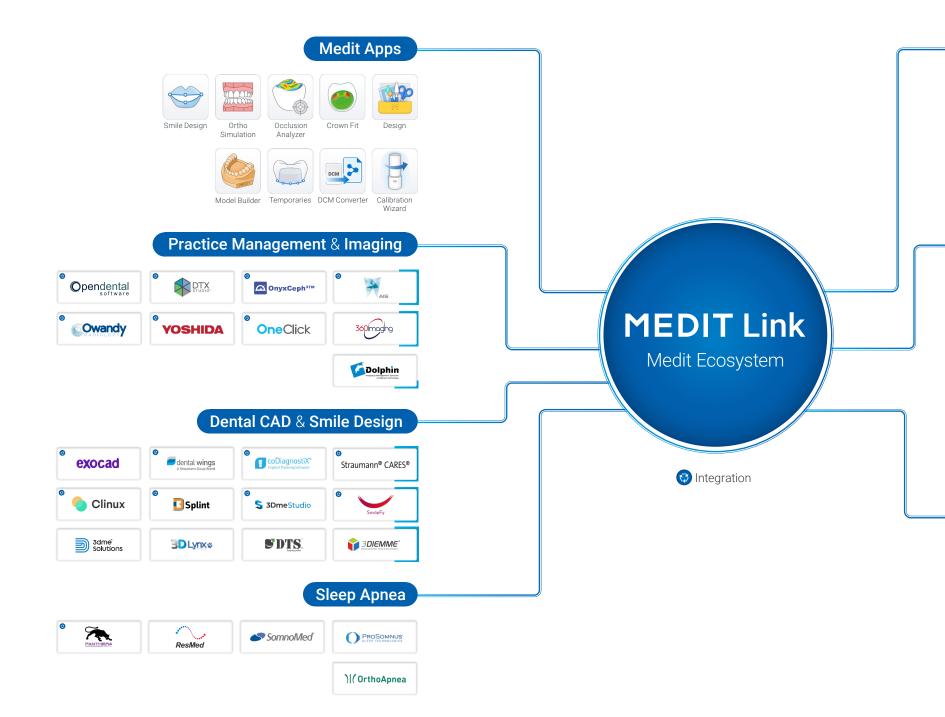

formlabs 🕅	8 ASIGA	SprintRay	<sup>®</sup> envision <b>tec</b>	ehitubox <sup>®</sup>	PYRAMIS	ACKURETTA	<sup>°</sup> Veltz	<sup>©</sup> ZENITH
Next Dent	HEYGEARS		ELEG	DEKEMA	🕖 phrozen	Phot <b>o</b> centric		W2P
<b>Denta</b> Fab	Sindoh 🏹		yndetect	SHERA Interest	carima 🞲	NOGASSAM		
Clear Alig	ners							
clearc <del>o</del> rrect	• reveal.	∞ ≫angelalign时代天使	● Smartee <sup>°</sup> 正册	KQ     Active Alignments	SP▲RK*		Science. Applied to Life."	<b>@ DynaFle</b>
		SMILE STYLER*	ACTIVE ALIGNERS	<b>Air</b> Smile®				
digitete ortho monufaktur		Clear Smile twister Mayer Clearly to M		Clast Smiletru Tixes & removable	crystalign	<b>DENTi</b> line <sup>®</sup>	Øexceed	
EXAKTUS	KLINE		휫 nutislab	orthoclear	ODL	ORGANICAL CAD CAM	<b>PRO</b> LIGNER	€e Profaligna
<b>₩₽</b> ₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽₽	sm;)ers	zeny∪m						
Others								
DDX	• easyrx	evident     Better by Design		MODIAW I	smile in a box.com	ABUTMENT	<b>_</b> \rchi <b>med</b> es <sup>,</sup>	<b>O</b> ALIENMILLIN
BlueSkyPlan.com	dicomL≜B <sup>°</sup>	DYNAMIC ABUTMENT « SCUTIONS			GuideMia Digital Dentitatry		C SmileWarp	Syntonia

# **Specifications**

Category		Medit <i>i</i> 700 wireless		
	Scanning Frame	Up to 70 FPS		
Coording Tools along	Imaging Technology	3D-in-motion video technology   3D full color streaming capture		
Scanning Technology	LED Source	LED		
	Anti-Fogging Technology	Adaptive anti-fogging		
Accuracy	Full-Arch	10.9µm ± 0.98		
l les defense	Dimensions	313 x 44 x 47.4 mm		
Handpiece	Total Weight	328 g (Includes battery, tip)		
	Tip Size	22.5 x 17.1 mm		
	Scan Area	15 x 13 mm		
Тір	Autoclavable	Autoclave 150 times 121°C 30 mins (Gravity Type) 134°C 4 mins (Pre-Vacuum Type) 135°C 10 mins (Gravity Type)		
	Reversible Tip	Yes		
	Remote Control Mode	Yes		
Special Mode	UV-C Disinfection	Yes		
Connectivity		Wireless USB 3.1 Gen1 (Wireless Hub)		



## About Us

# Our Medit Community





User Group

Help Center

Medit Academy

Medit specializes in 3D measurement and CAD/CAM solutions for dental clinics and labs, including intraoral scanners, based on its own patented state-of-the-art technology. We also develop platform solutions for digital dentistry, supporting collaborative workflows. Additionally, the Solutionix line provides 3D scanners and software to the industrial market. Our goal is to provide innovative technology and the highest quality products to ensure mutual growth for all partners.

Continuous feedback, suggestions of new treatment workflows, and active discussion. From hot debates to fun chats, the Medit community is always lively with our passionate customers. Colleagues from all over the world actively support and help one another, and sincerely welcome new colleagues into digital dentistry.

Scan the QR code to join!

Medit and its community of colleagues support and wish you the best of luck on your new journey into digital dentistry.

